BEBLISS

MyFirstScale

BABY DIGITAL SCALE

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■ DEAR USERS

Thanks for using our Smart Baby Scale. Please read and keep the manual in order to use the scale correctly.

This series of products can help you save and keep track of baby's weight, height and head girth on your smart device. With the APP, it can draw a growing curve versus WHO growing standard and , you can write a diary with photos to remind you of happy memories.

It can support multi-users.

■ FEATURES


1. Maximum capacity: 25kg/55lb; Division: 5g/0.01lb
2. Data comparison mark
3. Built-in measuring tape: 150cm / 60inch
4. Bluetooth connection, APP supported
5. App features: Measurement display, growing curve, diary and etc.

■ BUTTON FUNCTIONS

ZERO : On/Off ; ZERO

HOLD
UNIT : Short press to hold/unhold measurements,
long press to switch unit

■ BATTERY INSTALLATION

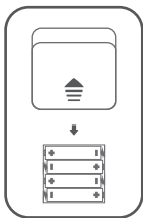Powered by 4 x 1.5V (AA/AAA) batteries .

First Use: Remove the insulation sheet before using (Note: If the battery is included)


Change battery: Take battery cover and install batteries as shown on the picture.

Note:

Please use the right battery.

Remove battery if the scale leaves unused for more than 3 months.


■ LCD ICONS INDICATION

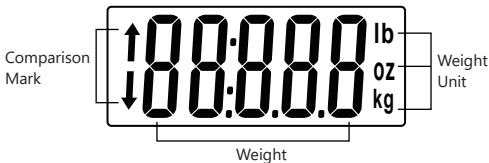■ WEIGHING


To ensure the accuracy of measurement, please use the scale on a flat and firm floor.

1. Press $\frac{\text{ON}}{\text{ZERO}}$ button to turn on and lay down the baby carefully on the scale.
2. LCD will show baby's weight and press $\frac{\text{HOLD}}{\text{UNIT}}$ to hold the measurement. The reading flashes twice for confirmation and will be stored in scale memory.
3. Comparison to previous data will be shown in difference value repeatedly. (Press $\frac{\text{HOLD}}{\text{UNIT}}$ to unhold or let it turn off automatically without any operation)

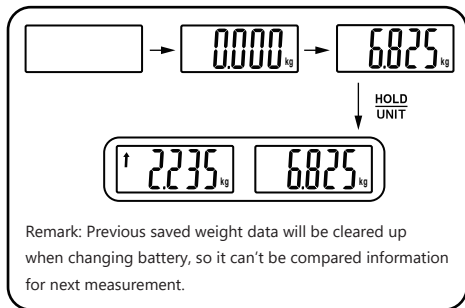■ MEASURING TAPE

The scale is with built - in measuring tape.
You can use the tape to measure the baby's head girth or height.

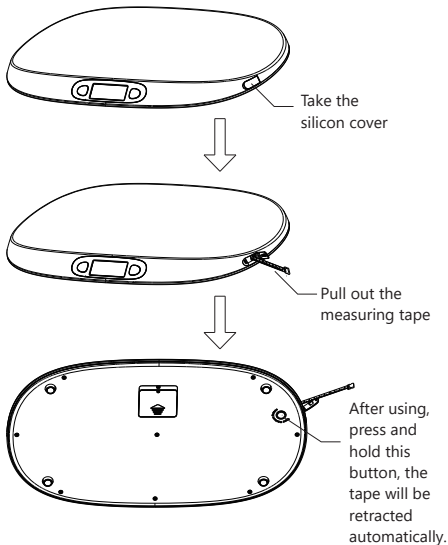
■ AUTO OFF AND MANUAL OFF

For longer battery life, the scale will turn off automatically after about 2 minute without any operation, or long press


"  " to turn it off after use.
ZERO

■ WARNING

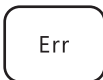OVERLOAD INDICATION

Scale overload the max capacity (Max capacity is indicated on the scale.) Please slowly hold up to avoid damage.

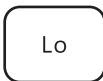LOW BATTERY INDICATION

Battery power is low.
Please replace new battery.

■ ADVICE FOR USE & CARE

1. Clean the scale with slightly damp cloth. DO NOT immerse the scale in water or use chemical/corrosive cleaner. Keep the scale in cool dry place.
2. All plastics parts should be cleaned immediately after contact with fats, spices, vinegar and strongly flavored/colored foods. (Avoid contact with acids citrus juice.)

3. Always use the scale on a solid, flat surface, keep the scale unmoved when weighing.
4. This product is precision instrument, please handle with care; Do not strike, shake or drop the scale.
5. If fail to turn on, please check whether battery(ies) are installed or battery power is low.
6. Keep away from high temperature, do not strip the scale arbitrarily.
7. This product can't be discarded as household garbage and should be sent to collection points for electrical and electronic equipment.
8. Designed for family use only, not for commercial use.

■ APP INSTRUCTION

1. Get the App

- a. Scan the QR code on the right side to obtain the App.
- b. At the "App Store" or "Google Play", search keyword "iMoni Health" to download.

APP: iMoni Health

2. Bluetooth Setting:

- a. Turn on the Bluetooth and Location on the smart device.
- b. Turn on the scale.

3. Create a new user

Run the App, finish the registration and personal profile.

4. Weigh and Measure

- a. Weigh the baby as direction, the weight will transfer to the APP.
- b. Measure the baby's head girth or height by using the built-in tape, and input the related measuring data to the APP.

5. Using the APP

- a. you can check the history of baby's weight.
- b. you can check the growing curve for baby within a period.
- c. you can write a diary with baby's photos.

For further information, please refer to the instruction in APP.

■ FAQ

The App can not synchronize the data from the scale

- Please ensure the connection between smart device and scale is on.
- Please ensure the Bluetooth function on smart device is on.
- Please ensure the scale is power on, and within the valid range of signal coverage.
- Check whether the battery power is running low, please replace with new batteries.
- Relaunch the App.
- Reboot the smart device and re-establish the connection with scale.

Bluetooth connection failed

- Please ensure only one device is setting up connection with scale.
- Please ensure your device is equipped with Bluetooth 4.0 or above.

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